



Water baths "Precisdig" with or without anti-evaporation folding lid

DIGITAL CONTROL AND DISPLAY OF TEMPERATURE AND TIME.

ADJUSTABLE TEMPERATURE RANGE FROM AMBIENT+5 °C TO 99.9 °C, STABILITY ± 0.1 °C, HOMOGENEITY ± 1 °C.

ADJUSTABLE TEMPERATURE RANGE FROM AMBIENT+5 °C TO 200 °C, STABILITY ± 0.5 °C, HOMOGENEITY ± 2.5 °C.

SAFETY: DIN 12877 STANDARD. THE UNIT WILL CUT OUT IF THERE IS INSUFFICIENT LIQUID OR IF THE BATH TEMPERATURE EXCEEDS THE SET TEMPERATURE. A MANUAL RESET ADJUSTABLE TEMPERATURE SAFETY THERMOSTAT IS PLACED AT THE BACK OF THE UNIT.

Programmable temperature and time. Interior and exterior made from Stainless steel




FEATURES

Double skinned tank of AISI 304 stainless steel with an external case of AISI 304 stainless steel.

Heating elements embedded in INCOLOY heat and corrosion resistant stainless steel.

Comes complete with drain tap.

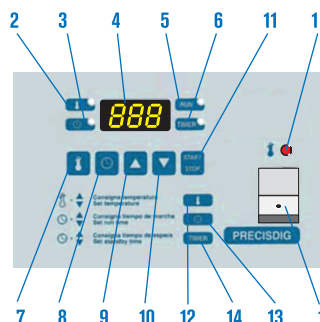

CONTROL PANEL

1. Illuminated mains switch.
2. Temperature mode indicator.
3. Time mode indicator.
4. Temperature and time indicator.
5. Start time indicator.
6. Waiting / delay time indicator.
7. Temperature set button.
8. Time set button.
9. Increase, push button.
10. Decrease, push button.
11. Stop and start time.
12. Read Temperature.
13. Set start up time: from 1 min. to 9 hours 59 min, or up to 99.9 hours once the set temperature has been reached.
14. Set waiting / delay time: from 1 to 24 hours.
15. Safety thermostat indicator lamp.

MAKROLON[®] FOLDING LID PROPERTIES

Due to its transparency, you can visualize the bath's inside. It improves temperature stability, and restricts evaporation and pollution

MODELS

Part No.	Temperature maximum °C	Capacity litres	Height / Width / Depth (usable) cm			Height / Width / Depth (exterior) cm			Power W	Weight Kg
6001195	99,9	5	15	30	14	20	37	27	600	6
6001236	200	5	15	30	14	20	37	27	1100	6
6001196	99,9	12	15	30	30	20	40	42	1100	8
6001237	200	12	15	30	30	20	40	42	1400	8
6001197	99,9	20	15	48	30	20	58	42	1500	10
6001238	200	20	15	48	30	20	58	42	2100	10
6001198	99,9	45	15	62	50	20	69	62	2400	18
6001239	200	45	15	62	50	20	69	62	3800	18

All baths with heating shelf and drainage. Faucet included.

MODELS WITH MAKROLON[®] FOLDING LID

Part No.	Temperature maximum °C	Capacity litres	Height / Width / Depth (usable) cm			Height / Width / Depth (exterior+lid) cm			Power W	Weight Kg
6001186	99,9	12	15	30	30	40	40	42	1100	8
6001187	99,9	20	15	48	30	40	58	42	1500	10
6001188	99,9	45	15	62	50	40	69	62	2400	18

All baths with heating shelf and drainage. Faucet included.

Marketed, Sales & Service by :

S.V. SCIENTIFIC Pvt. Ltd.

119, Sylvestra Arcade, 1st Floor, East Park Road, 11th Cross, Malleshwaram, BENGALURU - 560003. INDIA

Phone: 080-23567255 / 23462314 / 23446548

E Mail : info@svscientific.in / sales@svscientific.in

www.svscientific.com